VIBRATING NECK CUSHION 19958


CONTENTS

- 1 x Vibrating Neck Cushion
- 2 x AA 1.5V batteries

Please note: these batteries are basic batteries for demonstration purposes and may not last especially long – you may wish to replace them with better quality batteries)

BEFORE USE

1. Unpack – the cushion is supplied in a plastic bag, in a box.
2. Place the cushion on a clean surface so that the back of the cushion is nearest to you and the two ends are facing away from you.
3. Locate the zip at the back and open the zip.
4. Carefully remove the battery compartment and the batteries.
5. Open the battery compartment – pull back the tab and lift off the cover (there is no screw).

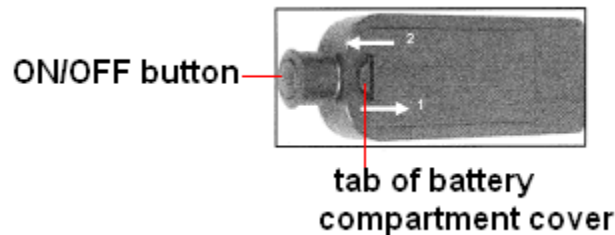6. Insert the two AA batteries, observing the polarity indicated.
7. Replace the battery compartment cover; first align the two 'pins' at the flat end and then gently push in the tab at the rounded end.
8. Press the circular button at the rounded end of the battery compartment to switch ON. A 'click' indicates when it has been pressed.
9. Place the battery compartment back in the cushion and close the zip.
10. Remember to switch off after use (press the circular button again).

CARE AND MAINTENANCE

Avoid sharp or pointed objects so as not to puncture the fabric.

Replace the batteries when operation becomes erratic.

Remove the batteries if the product will not be used for long periods.

Wipe the cushion with a clean sponge or cloth that has been slightly dampened with a mild cleaning product. Do not use a brush or anything abrasive. Do not use excessive amounts of water, benzene, window cleaning products, glass cleaning products, furniture polish or paint thinners. Pat dry gently with a dry cloth. Clean the battery holder with a cloth that has been slightly dampened with a non-abrasive, mild cleaning solution.

Do not wash

Do not bleach

Do not tumble dry

Do not iron

Do not dry clean

When this product is eventually exhausted, please dispose of it appropriately – not with domestic waste.

SAFETY

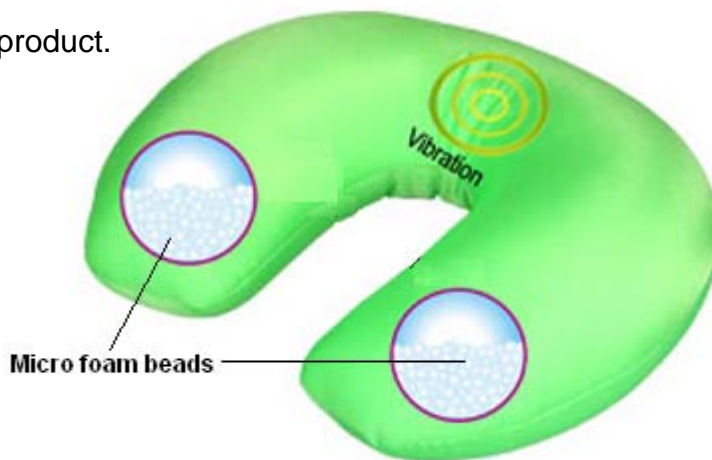
Always supervise the use of this product.

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